

iCOCOON



Relaxation Chair

HCP-ICO10000A

HCP-ICO10000D

HCP-ICO10000E

Safety precautions

Attached Document

Thank you for purchasing the iCOCOON.



This is a high-end, furniture-style relaxation chair that combines design and medication with the goal of promoting mind and body health.

Scan the QR code with your smart device to download the operating manual and the installation manuals.



Safety precautions







- Before using the product, read the following safety precautions so that you can use the product correctly.
- This reminder provides important information for the prevention of danger to the user of this product and surrounding people and property damage, as well as for the safe use of the product.
- Fully understand the following content (indication, pictogram), and then read the body text and observe the described matters.
- Descriptions of each chair component can be found on page 10-11.

	WARNING	This indicates a potentially hazardous situation, which, if not avoided, could result in death or serious personal injury.
	CAUTION	This indicates a potentially hazardous situation, which, if not avoided, could result in personal injury or property damage only.


 **WARNING**



- BURNS WILL RESULT FROM IMPROPER USE.
- READ INSTRUCTION BEFORE USING
- DO NOT WET - DO NOT USE PINS
- NEVER REMOVE COVER

Types and meanings of pictograms

 Prohibited	Indicates a prohibited act.	 Prohibited	Indicates that use in the bath, shower room, or other location is prohibited.
 Compulsory	Enforces and instructs an action.	 Disconnect plug	Instructs you to disconnect the power plug from the receptacle.
 No disassembly	Indicates that disassembly is prohibited.	 No wet hands	Indicates that handling with wet hand is prohibited.

Precautions on location of use

 **CAUTION**

 Prohibited	Do not use this product in a humid location such as a bathroom or sauna. → Failure to follow this instruction may cause an electric shock or failure.
 Compulsory	Install and use the product on a level floor surface. Installing the product on an unstable location may cause the main machine to fall, possibly leading to an accident.

After reading this reminder, keep it safe and available for future reference.

Precautions before use

WARNING

- The following persons must not use this product.

(1) People who are prohibited by a doctor to receive a massage

Examples: persons with symptoms such as thrombosis (embolism), severe aneurysm, acute varicosity, dermatitis, and skin infection (including inflammation of subcutaneous tissue)

(2) People who are suspected of having cervical sprain (whiplash)

- Do not use this product with a power supply other than the specified power supply.

→Failure to follow this instruction may cause a fire or electric shock.

Model number	Power supply
HCP-ICO10000A	AC120V ~
HCP-ICO10000D	AC220V · 240V ~
HCP-ICO10000E	AC110V ~

- Do not use this product by using a transformer.

→Failure to follow this instruction may cause a failure of this product or an electric shock.



Prohibited

- Do not use the product when the power cable or power plug is damaged or the power plug is loosely inserted into the receptacle.

→Failure to follow this instruction may cause an electric shock, short-circuit, or ignition.

- Before receiving a massage, be sure to check whether the cloth applied to the massaging portion is broken. Also check whether the cloth in other portions is broken. Even when you find a small breakage, immediately stop using the product, turn off the main power switch, disconnect the power plug from the receptacle, and then contact the store where you purchased the product. Never use the product by removing the cloth from the massaging portion or using broken cloth because doing so poses a great danger.

→Use of the product with cloth in a broken condition may cause an injury or electric shock.

- The following persons must consult a doctor before using the product.

(1) People who have heart trouble

(2) People who have sensory disturbance due to advanced peripheral circulation trouble that is attributable to diabetes or other

(3) People who have osteoporosis, broken spine, acute (painful) disease

(4) People who have a wound at the massaged portion

(5) People whose body temperature is 38°C or higher (duration of fever)

Example: period with strong acute inflammation symptom (malaise, chill, blood pressure fluctuation, and others), people in a debilitated condition

(6) People in an unstable period in the early stage of pregnancy, or people immediately after childbirth

(7) People with an in-vivo implanted medical electric device susceptible to electromagnetism, such as a pacemaker

(8) People who have a spine abnormality or spine curved to the right or left

(9) People who intend to receive a massage at a portion that was once treated by doctor or at a diseased portion that is under treatment

(10) People with a malignant tumor

(11) People who need bed rest

(12) People with acute symptomata such as visceral disease (gastric inflammation, intestinal inflammation, hepatic inflammation)

(13) People who are recognized to have lost temperature sensibility

(14) People under treatment by a doctor, particularly those who feel an abnormality in their body

(15) People who have deformed joints due to rheumatism, gout, or the like

(16) People with symptoms such as dizziness, lightheadedness, or tinnitus

(17) People with inflammation, such as a sprain or muscle strain

(18) People with tenosynovitis or a possibility of tenosynovitis






Compulsory

- Elderly people and persons with weak bones should seek medical advice before using the chair, even if they do not have a specific disease.

- If the use of the product does not have any effect, consult a doctor or specialist.

Continued from the previous page (Safety precautions)

Precautions before use	
 CAUTION	
 Prohibited	<ul style="list-style-type: none"> ● Do not use the product in a manner other than the specified manner. Do not use this product along with another medical device or an object such as an electric blanket. →Failure to follow this instruction may adversely affect your body or cause you to be injured.
	<ul style="list-style-type: none"> ● Do not carry the product with a person sitting on it. →Failure to follow this instruction may cause the main machine to fall, leading to an accident or injury.
	<ul style="list-style-type: none"> ● Do not sit on the chair with the foot rest raised. →Failure to follow this instruction may give a shock to the foot rest, possibly causing a failure. Failure to follow this instruction may cause the main machine to fall, leading to an accident or injury.
	<ul style="list-style-type: none"> ● Do not use the product by combining it with an item other than the supplied accessories. Do not use the supplied accessories by combining them with another medical device.
 Compulsory	<ul style="list-style-type: none"> ● Before inserting the power plug to a receptacle, confirm that the main power switch is turned off. →Failure to follow this instruction may cause a failure of this product or an electric shock.
	<ul style="list-style-type: none"> ● Positively insert the power plug and the connector all the way to the back. →Failure to follow this instruction may cause an electric shock, short-circuit, or ignition.
	<ul style="list-style-type: none"> ● Before using the product, turn on all the switches in order according to the separate installation manual, and confirm that the machine normally operates.
	<ul style="list-style-type: none"> ● During a massage, adjust the reclining angle according to your physical condition, and start massaging with weak stimulation. →If you start massaging with strong stimulation, you may be injured. For elderly people or people with weak bones in particular, massaging must be performed with sufficient care.
	<ul style="list-style-type: none"> ● For your safety, check the positions of the mechanisms. If they are not at their housing positions, turn on the main power switch, confirm that the machine operates normally, and then sit on the chair.
	<ul style="list-style-type: none"> ● When sitting on the chair, do not step, stand, or sit on the foot rest. →Otherwise, an accident or injury due to falling or failure may be caused.
	<ul style="list-style-type: none"> ● When you are going to use the product again after a long period of disuse, carefully read this instruction manual again, confirm that the machine normally operates, and then use the product.
	<ul style="list-style-type: none"> ● Before using the product, be sure to completely raise the back rest and completely lower the foot rest.

Precautions before use

WARNING



Prohibited

- Do not wind the power cable around the main unit, or step on it with the leg of the main machine. Do not damage, break, process, forcefully bend, pull, or twist the power cable.
→Failure to follow these instructions may break the cable, leading to a fire or electric shock
- Do not spill water or other liquid on the main unit, main unit remote control or monitor.
→Failure to follow this instruction may cause an electric shock, short-circuit, or ignition
- Make sure that children, people who cannot express their will, or physically-handicapped people do not use the product. Do not let a child or a pet play on the main machine or around the main machine (behind the back rest, beneath the foot rest, etc.), climb on the main machine, or sit on it. Only in the meditation program and the child program, a child of 10 years or older may use the product. In this case, a guardian must always attend the child.
→Failure to follow this instruction may leading to an accident or injury.
- This appliance can be used by children aged from 10 years and above and persons with reduced physical, sensory or mental capability or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved..
→May result in accident or injury.
- Children shall not play with the appliance.
→May result in accident or injury.
- Do not use the product by letting two or more persons sit on the chair. In particular, do not use the product while holding a child (baby, etc.).
→Failure to follow this instruction may leading to an accident or injury.
- Do not use the product with the cloth in a broken condition. Disconnect the power plug from the receptacle.
→Failure to follow this instruction may case an accident, injury, or electric shock. .
- When you are holding the foot rest with your hands, do not suddenly release your hands.



No wet hands

- Do not sit on the chair with a wet body or operate the product with a wet hand.
→Failure to follow this instruction may cause an electric shock or failure.



Compulsory

- Before you start a massage course or when reclining the back rest and returning to the original posture, be sure to confirm that there are no persons (children in particular) and no pets around you. Do not have your body or an object caught at the back of the back rest or beneath the foot rest.
→Failure to follow this instruction may lead to an accident or injury.
- When your body or an object gets caught, press the reclining button to return the back rest or foot rest to a safe position.

CAUTION






Prohibited

- Do not use the product with bare skin.
- Do not use the product for purposes other than massage and sitting.
→Failure to follow this instruction may lead to an accident or injury.
- Do not use the product on the head, chest, stomach, and joints (elbows, knees, and others).
→Failure to follow this instruction may cause you to feel sick or get injured.
- Do not put your hand or foot on the path where the massage roller moves up and down, nor on movable parts such as the reclining back rest and mechanisms. Do not use the product when wearing hard objects such as hair accessories, piercings, earrings, a necklace, a wrist watch, bracelets, finger rings, and accessories.
→Failure to follow these instructions may cause you to get injured.
- Do not use the arm rests on other body parts than arms.
→Failure to follow these instructions may cause you to get injured.

Continued from the previous page (Safety precautions)

CAUTION

 Prohibited	<ul style="list-style-type: none"> ● Do not use the foot rest on body parts other than the feet. →Failure to follow these instructions may cause you to get injured. ● Do not operate the massage chair with no one present. →Failure to follow this instruction may lead to an accident or injury
 Compulsory	<ul style="list-style-type: none"> ● In an emergency situation or when you have sensed any abnormality, press the [Stop] button on the tablet to raise the reclined back rest ● When you feel a strong pain in your body or sense any abnormality during use of the product, immediately stop using it and consult a doctor. ● If the massage stimulation is too strong for you, weaken the stimulation by following the instruction manual. ● If the stimulation is still strong, immediately stop using the product and consult the store where you purchased the product. →Continued use of the product with a strong stimulation may cause you to get injured. ● During a massage, be careful not to fall asleep. →Failure to follow these instructions may cause you to get injured. ● The heater is located under the foot and back of the chair. Person who cannot feel heat should be careful when using the heater. ● The appliance has heated surface. Persons insensitive to heat must be careful when using the appliance. ● Do not receive massage for over 30 minutes a day. Continued use of the product on the same portion must be for a maximum of 5 minutes. → Continued use of the product for long hours provides stimulation that is greater than necessary, which may cause an adverse effect and/or injury. Receive massage once a day as a guide. ● If your arm has come out of the arm rest during massaging of the arm, do not forcefully put your arm back in the arm rest. →Failure to follow this instruction may cause failure.。 ● If your foot has come out of the foot rest during massaging of the foot, do not forcefully put your foot back in the foot rest. →Failure to follow this instruction may cause failure. ● When you receive a massage at the area around the neck, pay attention to the position of the massage ball so as to prevent massaging of the throat and the front portion of the neck. ● Do not turn up the volume to maximum level from the beginning. Make sure not to turn up the volume of the external music system connected to the speaker audio input from the beginning. Doing so may produce an unexpected, loud sound at the time of turning on the speaker, causing damage to your hearing.Failure to follow this instruction may cause failure.
 Disconnect plug	<ul style="list-style-type: none"> ● When a blackout has occurred or there is a possibility of blackout, immediately stop using the product, turn off the main power switch, and disconnect the power plug from the receptacle. →When the power is restored, there is a risk of an accident or injury. ● When there is a possibility of lightning strike, immediately stop using the product, turn off the main power switch, and disconnect the power plug from the receptacle. ● If the product does not operate or you have sensed any abnormality, for the purpose of preventing accident, immediately stop using the product, turn off the main power switch, disconnect the power plug from the receptacle, and request for the store where you purchased the product to inspect and repair it. →Failure to follow this instruction may cause an electric shock, electric leakage, short circuit, or other defect, possibly leading to fire.

Precautions after use or during disuse of product

CAUTION



Prohibited

- Do not use the power cable on a device other than this product.
- After using this product, be sure to completely raise the back rest and completely lower the foot rest.
- After using the product, wipe off any dirt, and store the product in a location subject to little humidity.
- Do not put an object on the product.
- Store the cables untangled.



Disconnect plug

- After use of the product, be sure to turn off the main power switch, and disconnect the power plug from the receptacle.
→ Failure to follow this instruction may cause an accident or injury due to mischief made by a child or pet.
- When disconnecting the power plug from the receptacle, be sure to hold the power plug, not the power cable.
→ Failure to follow this instruction may cause an electric shock or short circuit, possibly leading to fire.
- When taking care of the product, be sure to disconnect the power plug from the receptacle. Do not insert or disconnect the power plug with wet hands.
- When you are not using the product, be sure to disconnect the power plug from the receptacle.
→ Failure to follow this instruction may cause dust and humidity to deteriorate the insulation, possibly leading to an electric shock, electric leakage, or fire.

After use or when not in use

WARNING



Compulsory

- Cleaning and user maintenance shall not be made by children without supervision.
→ May result in accident or injury.

Other precautions

WARNING



Prohibited

- Do not stand on the main machine. Do not give a strong shock to the main machine from above by, for example, sitting on a portion other than the seat.
→ Otherwise, an accident or injury due to falling or failure may be caused.
- When the power cable has been broken, immediately stop using the product because use of the product in this condition poses a danger.
→ Request for the store where you purchased the product to replace the power cable.



No disassembly


- When the machine has failed, do not repair it by yourself.
→ Contact the store where you purchased the product.
- Do not modify the product.
→ Failure to follow this instruction may cause a fire, electric shock, or injury.






Request

When disposing of the product or its accessories, contact the nearest municipality office for the correct disposal.

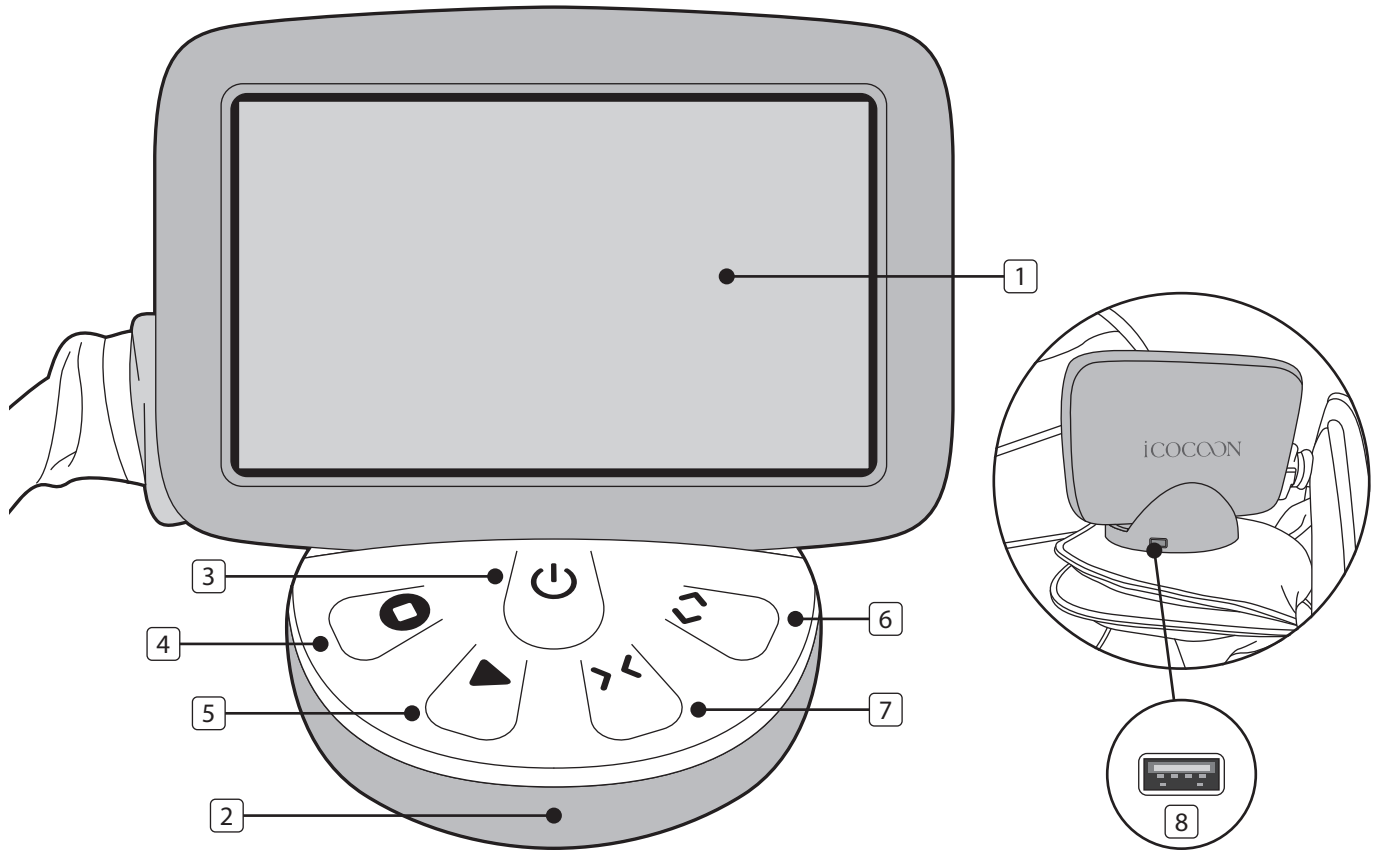Names and functions of machine components


Monitor / Main Unit Remote Control



①	Monitor (Accessory)	Operate by tapping the screen.
②	Main Unit Remote Control	You can operate the main unit and also can install the monitor.
③	 ON/OFF Button	Turns the main unit power "ON/OFF." When the remote control power is "ON," it lights up green. When it is "OFF," the light turns off. When the power turns off, the reclining seat back will rise and massage rollers will return to the storage position before stopping. The reclining seat back automatically rises when pressing the button during the stage of emergency stop.
④	 Emergency Stop Button	Press this button in case of an emergency or if you sense any abnormality during the operation. The operation will stop immediately. When an emergency stop is being made, the button flashes in red.
⑤	 FULL BODY	This is a recommended program that provides a full-body massage that helps with recovery from fatigue.
⑥	 Open Button	The side panels will open. They will continue opening while holding down the button.
⑦	 Close Button	The side panels will close. They will continue closing while holding down the button.
⑧	USB Port	This is a communication port for customer support. Please do not use it.

Continued (Name and functions of machine components)

Main Unit (Front side)

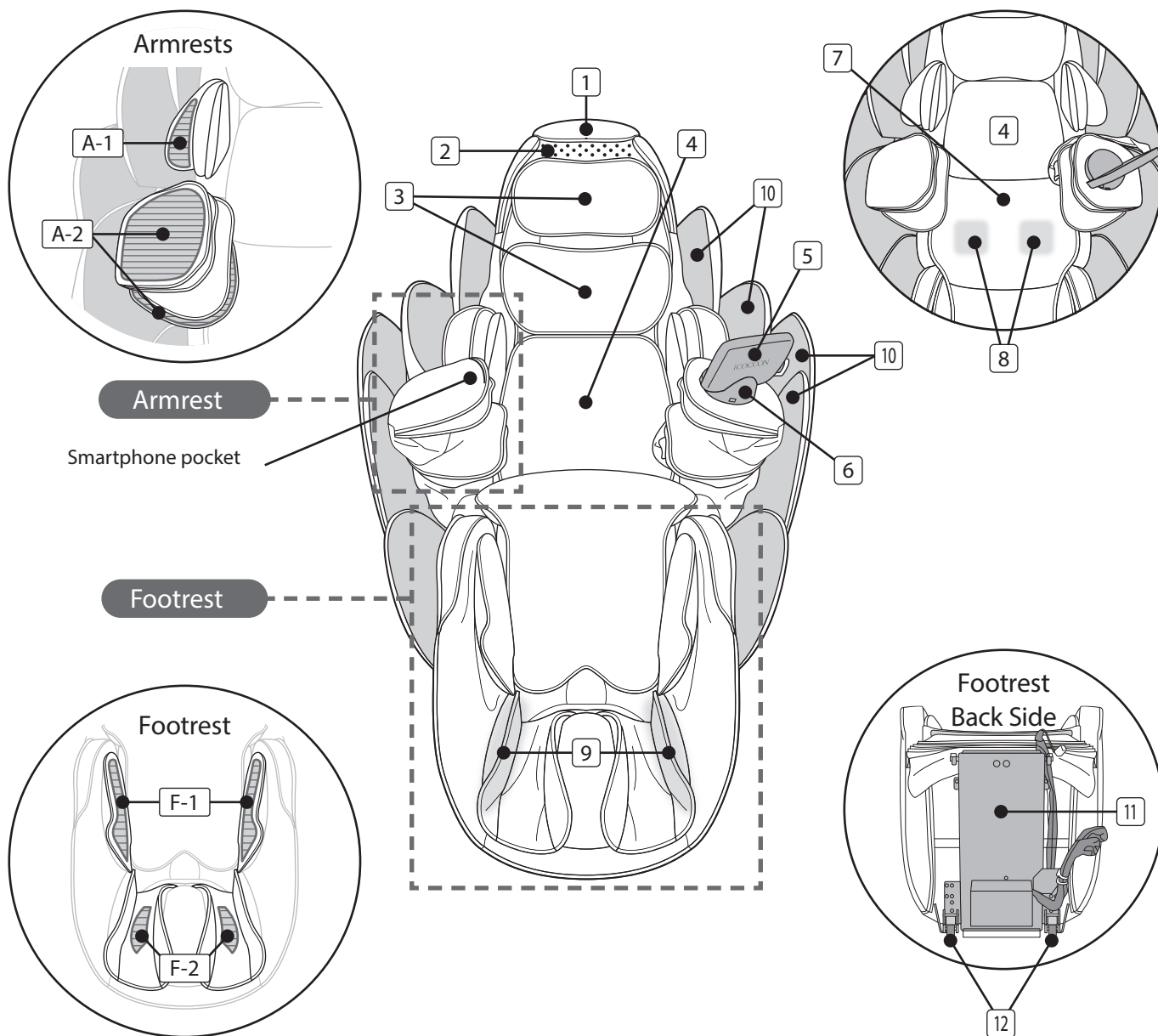

1	Top Panel
2	Speakers
3	Pillow
4	Back Cushion

5	Monitor (with Monitor Cable)
6	Main Unit Remote Control
7	Seat Cushion
8	Seat Heater

9	Foot Heaters(right and left)
10	Side Panels(right and left)
11	Foot Slide
12	Casters (right and left)


Armrest

A-1	Shoulder Mechanism(right and left)
A-2	Arm Mechanism (right and left)

Footrest

F-1	Calf Mechanism (right and left)
F-2	Foot Sole Mechanism (right and left)

Main Unit (Back side)

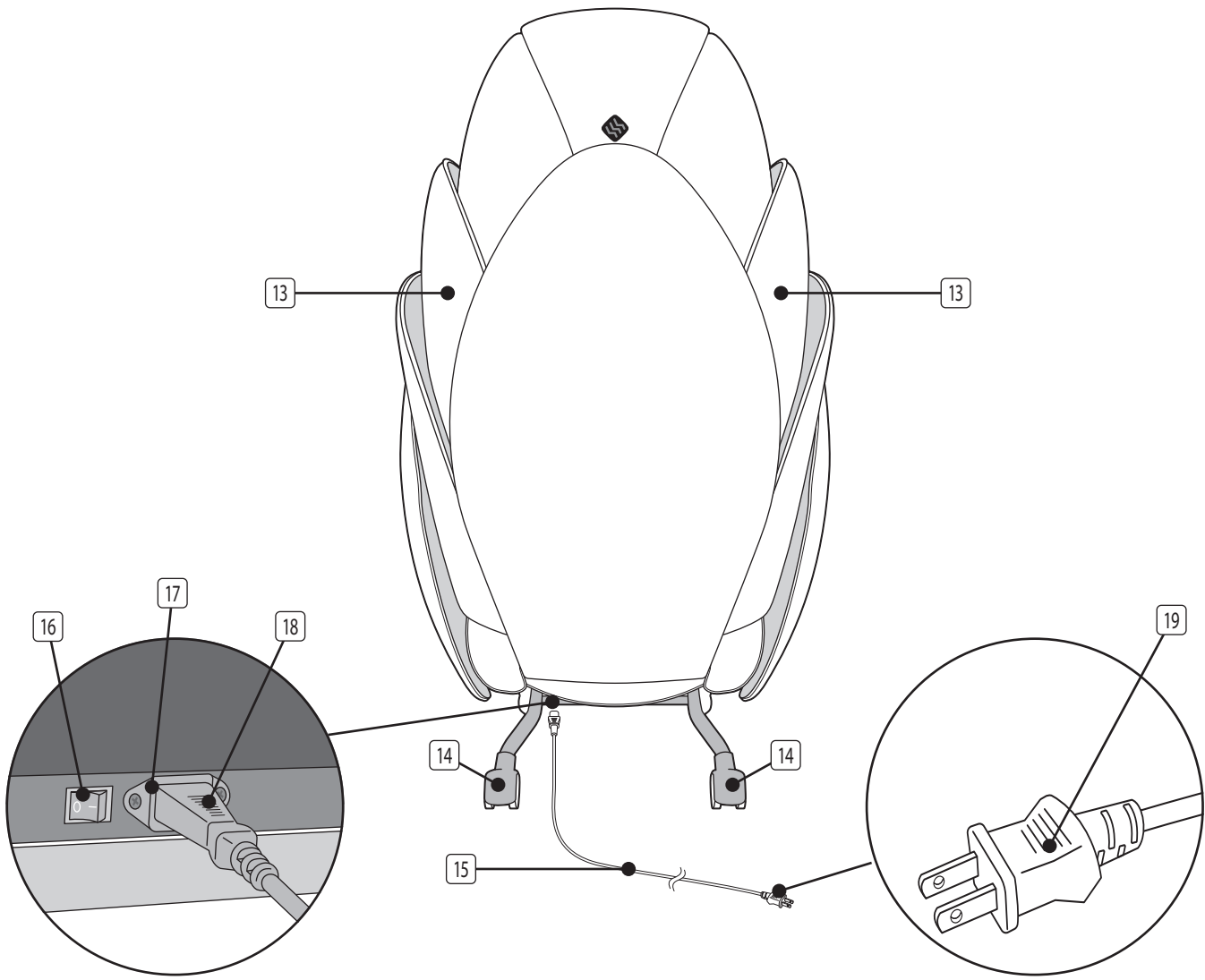

13	Side Panel (right and left)
14	Caster (right and left)
15	Power Cord
16	Main Power Switch


17	Inlet
18	Connector
19	Power Plug

Air Cells

K-1	Pillow Air Cell
K-2	Lower Back Right Air Cell
K-3	Lower Back Left Air Cell
K-4	Seat Right Air Cell

K-5	Seat Left Air Cell
K-6	Seat Front Air Cell
K-7	Back Cushion Air Cell
K-8	Seat Cushion Air Cell


Air Cells

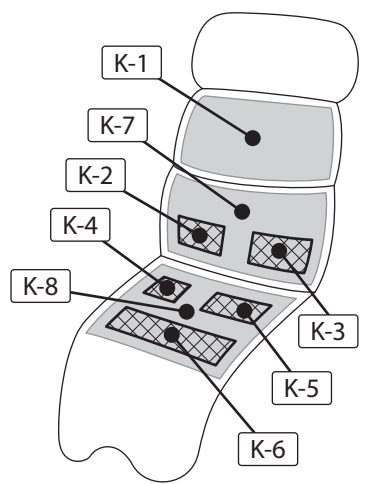

Before Starting Massage

Check around the chair

Make sure that no one, pets or objects are around or behind the chair or in front of the leg section.

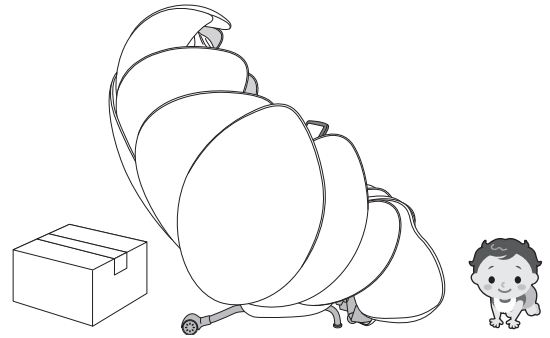

Turn on the power

Warning

- Use the product only at the specified power source.
→Using the wrong power source may result in a fire or electric shock.
- Do not use the product with a transformer. This may result in malfunction or electric shock.
- Do not use the product if the power cord or power plug is damaged or its insertion into a wall outlet is loose.
→This may result in an electric shock, short circuit, or fire.

Caution

- When using the product after raising the ambient temperature in a very cold room, do not increase the room temperature suddenly but gradually.
→Condensation may occur on the operating parts, and performance may not be fully demonstrated.


Main Power Connection

- 1 Plug the power cord into the cord connector located at the back of the chair first. Then plug the other end into the wall outlet.



Set the main POWER switch located at the back of the chair to ON.

A text message should appear on the remote control display screen.

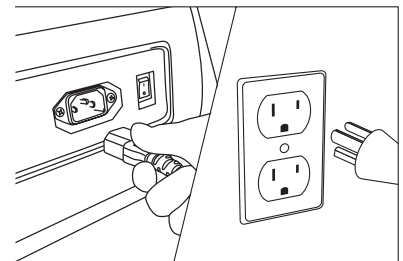

* If the main POWER switch is turned on while the rollers are stopped in the middle of the movement, they will automatically return to their storage position (upper part of the reclining seat) then stop. The seat will also move for adjustment of the initial position then stop.

Warning

- Do not use with any type of transformer device.
→May result in damage to the product or electric shock.


Caution

- Before plugging the power cord into the wall outlet, make sure that the main POWER switch is set to OFF.
- Make sure that the plugs on the power cord are inserted fully into the cord connector and the wall outlet.
- Before using the chair for the first time, turn on all the switches in order as indicated in this manual to check that the chair is operating properly.


Continued (Before Starting Massage)

2 Turn **ON** the Main Power Switch of the main unit.

- The main power switch is located on the back side of the main unit.

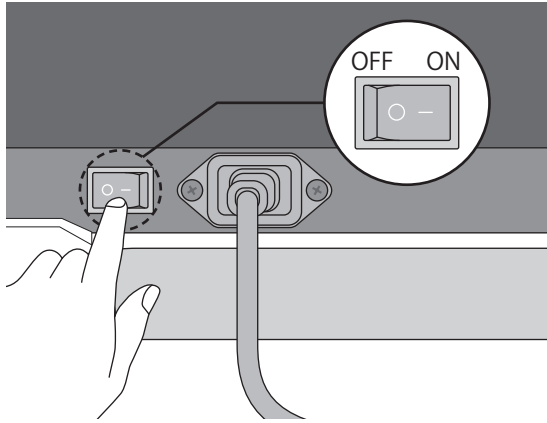

3 The startup screen appears on the monitor screen. During preparation, the LED on the **ON/OFF** button flashes in green.

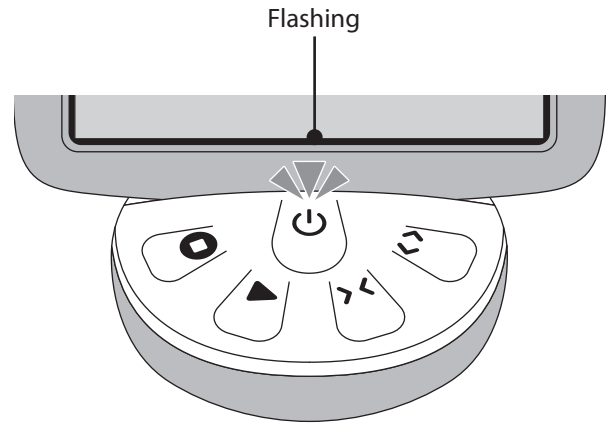

Check the main unit

1 Make sure to check that the fabric of the main unit is not torn.

Warning

- Before use, make sure to check that the fabric of the main unit is not torn. Also, make sure to check that the fabric of other parts is not torn.
- However small a break, immediately stop using the product, disconnect the power plug, and contact your local distributor. Never remove the fabric of the main unit or use it in a torn state as it is very dangerous.

2 Make sure that no foreign objects are caught in the main unit.

Caution

- Make sure to check that no foreign objects are caught in the main unit before use.
→May cause an accident or injury.

3 Make sure to check that each mechanism is in its stowed position.

Mechanism stowed positions

- Back mechanism ... Top of the backrest
- Sole mechanism ... Open state

Caution

- Make sure that each mechanism is in the stowed position before sitting down.
→May cause an accident or injury.

4 Make sure to check that the backrest is fully raised and the footrest is fully lowered.

- * When the power is on, the side panels will open, and air will be supplied to the pillow, back cushion, and lumbar cushion.

Caution


- When the footrest is not fully lowered, do not sit down, stand, or place an object on the footrest.
→May cause an accident or injury.

How to sit down on the chair

Normally, each mechanism is in a stowed state.

Back mechanism ... Bottom of the backrest; Sole mechanism ... Open state

For safety, check each mechanism position before sitting down, and if it is not in a stowed state, do not sit down immediately but turn on the main power switch to return each mechanism to the stowed state, and then slowly sit down in the center of the chair.

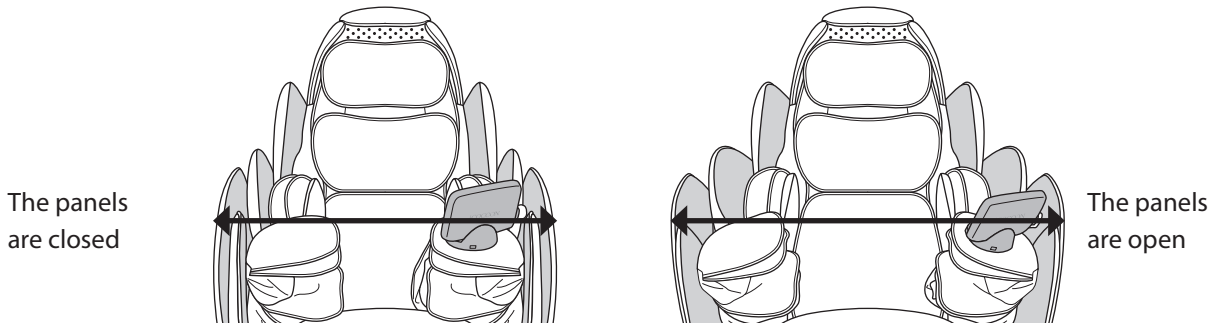


Caution

- Do not sit in the chair while the footrest is raised. An impact may be applied to the footrest causing a malfunction. Also, an accident or injury may result in the case of the main unit overturning.
- Before use, fully raise the backrest.

Precautions for getting in/out of the chair

- 1 Make sure to check that the footrest is fully lowered, and then put your feet into the sole unit.
- 2 Get in/out of the chair by placing your hands on the upper part of the footrest.

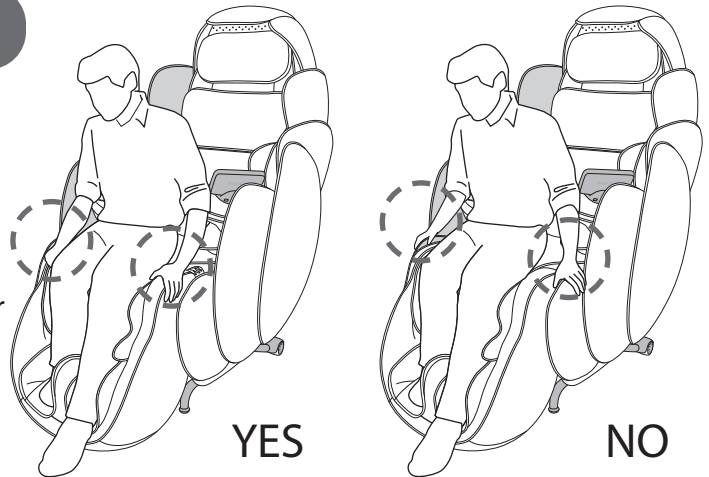

Warning

- Do not subject the side panels, main unit remote control or armrests to strong impacts or lean against them. →Overturning them may cause an accident, injury, or a malfunction.
- If the arm mechanism and/or footrest is in operation, wait for the operation to stop, and then get in/out of the chair.
- Be sure to sit in the chair with the side panels open.

Caution

- Do not sit down on the footrest while it is raised.
- Never get in/out while each mechanism is moving because it is very dangerous. →May result in injury or a malfunction.
- We recommend use of this product for persons weighing 100kg (220 lbs.) or less. When a person weighing 100kg (220 lbs.) sits in the chair, a load of approximately 1.9kg (4.2 lbs.) is applied to the floor per 1cm² (0.16 in²).
- Before use, make sure to check that the fabric is not torn.
- The withstand load of the armrests is approximately 30kg (66 lbs.).
- The tensile load of the armrests is approximately 30kg (66 lbs.).
- The withstand load of the footrest is approximately 30kg (66 lbs.).

Continued (Before Starting Massage)

Arm positions

Put your arms in between the armrest.

- Do not subject the armrest to strong impacts or lean on it. It may cause deformation or damage.

Correct arm positions



Place your hands palm-down.

Warning

- Do not place your arms outside the elbow rests of the armrests.
- The arm mechanisms slide. Do not put your hands or fingers in the spaces between the arm mechanisms and armrests.
→May result in injury.


Caution

- Do not use the armrests for anything other than your arms.
→May result in injury.
- If one of your arms comes out of the armrest while your arms are being massaged, do not try to force it back in. Undue force may cause a malfunction.
- Do not massage your elbows.
→May result in injury.

Leg positions

1 Put your legs in between the calf mechanisms and footrest as shown in the illustration.

2 Sit back in the seat. It does not matter if your feet are not touching the bottom of the footrest.

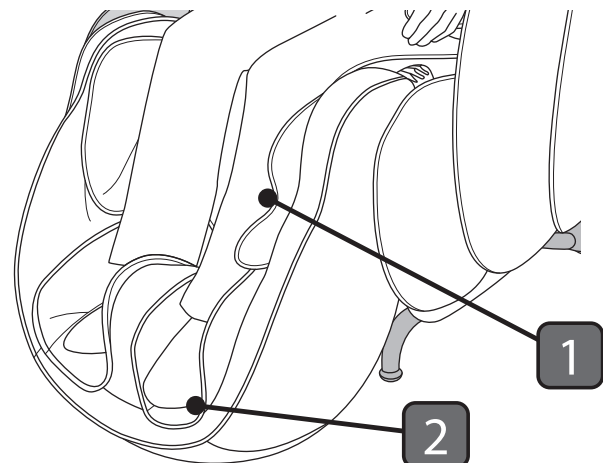


Automatic Foot Slide Adjustment Function

The Automatic Foot Slide Adjustment Function extends and retracts the footrest by responding to a sensor in the toe area of the footrest. For sensors to react properly, sit down with your soles firmly pressed against the footrest.

* If your toes do not reach the footrest, make adjustments to an angle of the backrest until your legs fit the footrest. The foot slide will be automatically adjusted when starting the massage program.

Extension and retraction of the footrest

When raised to a certain height, the footrest can be extended/retracted.

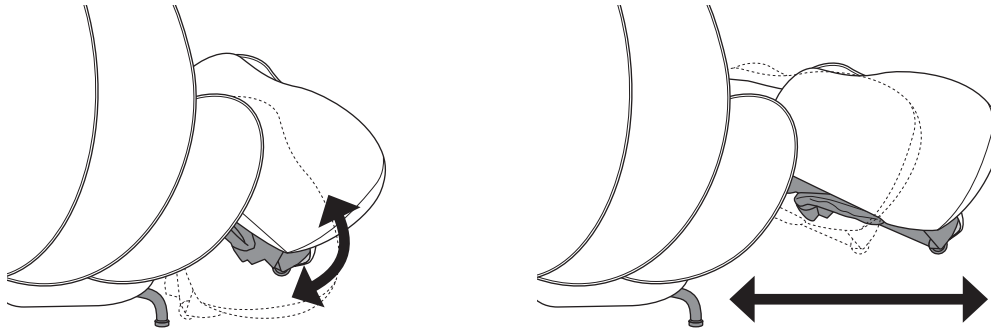

*The extending/retracting length is approximately 27cm (11 inch).

Warning


- Do not sit on the footrest while it is being raised.
- When lowering the footrest, make sure to check there is no one (child in particular) or any pet around the chair. Also, be careful not to get your body or other objects caught under the footrest.
- Do not suddenly release your hand when lifting the footrest.
- When the foot slide function is in motion, be careful not to get your feet or other objects caught in between the foot slide and the footrest.

Caution

- Do not use the footrest for parts other than your legs.
→May result in injury.
- If your one of your legs comes out of the footrest while your legs are being massaged, do not try to force it back in.
Undue force may cause a malfunction.
- Do not massage your knees. Do not put your legs in other directions.
→May result in injury.

Monitor Position

Remove the monitor located between the back cushions and attach it to the main unit remote control. The monitor and the main unit remote control are attached together with magnets.

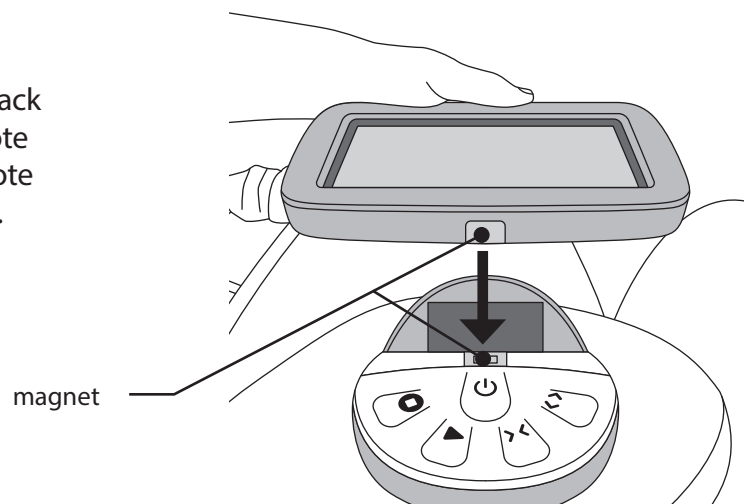


Caution

- Do not use if you have an electronic medical device such as a pacemaker.
- Do not bring any magnetic storage media, such as magnetic cards, near the chair. It may result in destruction of data. Additionally, do not bring PCs, electronic wristwatches, or other precision electronic devices close to the chair, as this may cause them to malfunction.
- Do not move the angle adjustment part of the monitor cable beyond its movable range or try to forcibly pull or bend it.
→ The harness inside the monitor cable may break and cause damage.
- Do not apply any load to the monitor or the main unit remote control, such as leaning against or giving a strong impact to them.
→ This may result in an accident, injury or malfunction.


How to Use the Monitor (Before sitting in the iCOCOON)

How to activate the monitor

1 Turn **ON** the main power switch located at the back left of the iCOCOON main unit.

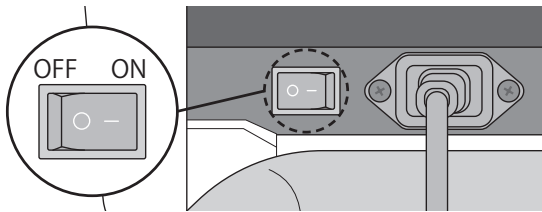

The startup screen appears on the monitor screen. During preparation, the LED on the **ON/OFF** button flashes in green.

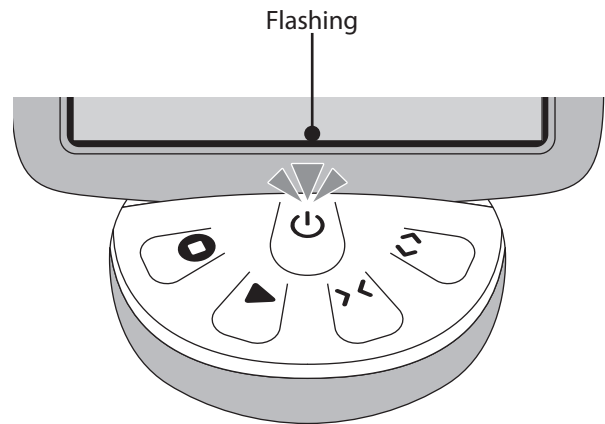


2 When the standby screen appears, press the **ON/OFF** button of the main unit remote control.

- After the elapse of a certain period of time with no operation being performed, the screen turns black. press the **ON/OFF** button to turn it on again.

* For details, refer to the FAQ.

3 The power is turned "ON" and the iCOCOON panels open. When the power is ON, the LED on the **ON/OFF** button lights up green. After that, air is automatically supplied to the pillow, back cushion and seat cushion.

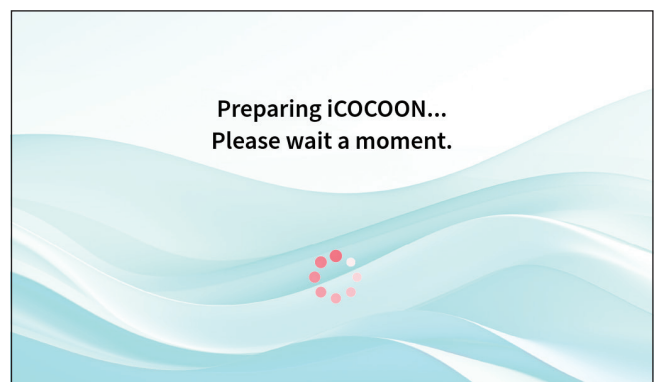

4 The User Selection screen will appear after air preparation is completed. Please sit in the iCOCOON.

